# The solution to manage low volume immunoassay testing.





### Reliability


- Utilises the proven AIA-PACK concept
- · Positive identification of samples and reagents
- Fibrin sample clot, air bubble detection and sample level sensing

# **Operational Simplicity**

- Start up in less than 10 minutes
- · Limited maintenance
- No need for test programming. Simply place the samples and AIA-PACKs and press START

# **Flexibility**

- Same reagents used across the entire AIA franchise
- Its combined features make it an ideal instrument

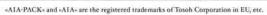





EXCLUSIVE DISTRIBUTOR FOR MONTENEGRO

Mrkosnica 18 - 81400 Niksic - MONTENEGRO Tel/fax: +382 (0) 40 214 500 - info@super-lab.me www.super-lab.me




Transportstraat 4 - 3980 Tessenderlo - BELGIUM Tel: +32 (0)13 66 88 30 - Fax: +32 (0)13 66 47 49 www.tosohbioscience.eu



### TOSOH CORPORATION - BIOSCIENCE DIVISION

Shiba-Koen First Bldg. 3-8-2, Shiba, Minato-ku, Tokyo 105-8623 - JAPAN Tel: +81-3-5427-5181 - Fax: +81-3-5427-5220 www.tosoh.com

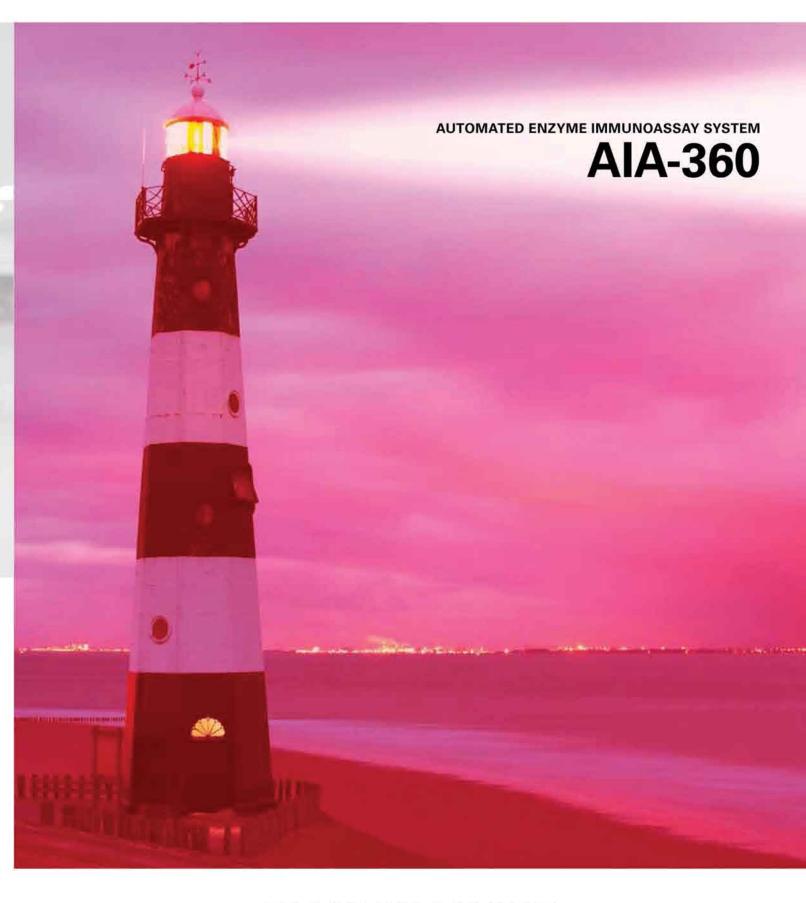

### **TOSOH BIOSCIENCE**


## The AIA-PACK. A unique concept.

- Lyophilised reagents allow long AIA-PACK shelf-life\*
- 90 days calibration stability
- Most of the analytes require a 10 minute incubation time, delivering the first result within 20 minutes
- Low risk for contamination or carry-over between samples and reagents:
  - · I AIA-PACK = I test
  - The complete assay (incubation, wash and detection) is carried out in the AIA-PACK
- The single cup per test combined with low volume consumption of reagents reflects Tosoh's environmental awareness

<sup>\* 12-</sup>months from date of manufacture



A ready-to-use "all-in-one" test cup

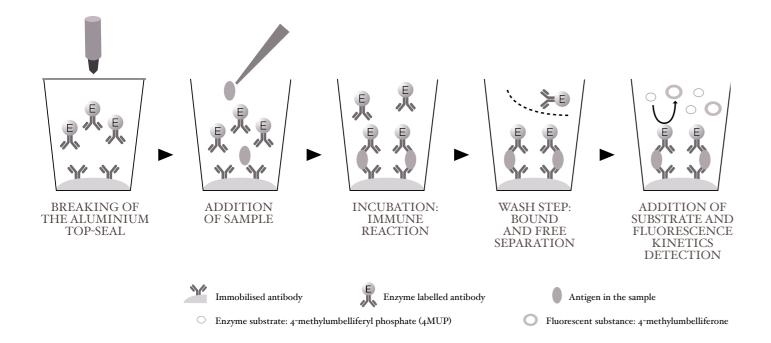

# A real throughput of 36 tests/hour with continous loading capacity.

### Friendly operator interface

- The AIA-360's icon driven touch screen is intuitive and allows for easy navigation
- Non bar-coded specimens can be manually programmed
- Up to 800 assay results can be saved on the AIA-360
- Traceability of all samples and data
- Onboard full QC package including Levey-Jennings charts\*
- 2 calibration curves simultaneously active per analyte
- Recalculation of results in case of recalibration
- Fully customisable patient report\*
- Patient results printed in real time for fast and easy review
- Interfacing allows connection to Laboratory Information Systems (LIS)\*

\*In combination with TOSOH EUROPE N.V. PIANO software



# #IIA-360

### Advanced sample carousel

- Up to 25 patient samples and AIA-PACK tests in a single run
- Primary tubes as well as sample cups can be loaded onto the system
- Random access: continuous loading capacity of samples
- Continuous loading capacity of reagents

### Common reagents and waste area

- AIA-PACK principle allows most of the waste to be treated as solid waste with a significant reduction of disposal costs
- Sophisticated management system based on level sensing and alarms for reagents and waste

Dimensions: W 400 mm, D 400 mm, H 520 mm, Weight: 29 kg